





#### **SURGICAL DIODE LASER**

CONTROL THE BLEEDING



CONTROL THE OUTCOME **(** 





## VET REVIEWS







"Most user friendly laser on the market, much faster sarcoid removal compared to competitor machine we used to hire!"

Jess Watson - Hampden Equine



"MED Equus' diode laser machine has been revolutionary for our large oncology caseload. The machine is easy to set up and the ability to change settings according to tissue type ensures exceptional control for sarcoid and melanoma removal and for lack of bleeding and haemostasis.

KENT EQUINE HOSPITAL

Katie Snalune - Owner and Head of Kent Equine Hospital

# VET REVIEWS





"We purchased the Leonardo laser because we trust MED Equus to provide high quality equipment and good after sales support."

James McIntosh - Southmoor Equine



"I like the user interface with ease of adjustable settings and portable. Being able to adjust the setting to vary the proportion of each wavelength is useful. For delicate areas we use less of the 980nm wavelength to reduce thermal damage."











# THE COMBINED WAVELENGTH WAS REALLY EFFECTIVE DURING OUR FIRST 3 MONTHS OF SURGERY

Leading charity completed 16 surgeries in 3 months and only had 1 case that bled



OLD V NEW











#### **Leonardo - Technical Specification**

Wavelength	980nm and 1470nm	Pulse duration	0.01 - 60sec / 0.01 -60 sec
Power	Max. 45WATTS	Power supply	110 - 24-VAC, 50/60Hz, 450VA
Fibre diameter	≥ 360 µm	Weight	8.5kg
Aiming beam	532nm and 635nm, green 1mW, red 4mW, user controlled intensity	Dimensions	28Hcm x 37Wcm x 9Dcm

**Treatment mode** 

CW, Continuous and pulse

## Lemnardo DUAL 45

- 980nm and 1470nm wavelength combined for optimum tissue interaction
- Increased power output to 45W
  Intuitive, touch screen operation
- Transparent safety goggles
- Green aiming beam instead of red for ease of use



45 WATT - German engineered



Dual wavelength - a quicker and cleaner cut



Adjustable modes for airway and skin applications



Average product life span over 10 years



Unique haemostasis setting for vessel sealing



## Lemordo DUAL 45

#### PRE-PROGRAMMED SETTINGS

#### - SARCOID 1

PRECISE SETTING FOR REMOVAL OF SARCOIDS FROM EARS, NEAR EYES AND VASCULAR AREAS

#### - SARCOID 2

PRECISE SETTING FOR DE-BULKING SARCOIDS/MELANOMAS

#### - HAEMOSTASIS

ENHANCED SETTING TO STOP BLEEDING OF A VESSEL REDUCING THE NEED TO USE A LIGATURE

#### - SCARIFY

OPTIONAL REDUCED POWER SETTING TO SCARIFY THE AREA WHERE THE SARCOID HAS BEEN REMOVED

#### **AIRWAY RECOMMENDATIONS**

DEPENDING ON THE CLINICAL AIRWAY PROCEDURE YOU ARE PREFORMING THE LEONARDO IS THE PERFECT SURGICAL TOOL.

THE ABILITY TO HAVE THE COMBINED WAVELENGTH IS CRUCIAL WHETHER YOU ARE CARRYING OUT A VENTRICULAR CORDECTOMY, RESECTING ARYEPIGLOTTIC FOLDS OR REMOVING AN EPIGLOTTIC CYST.

THE COMBINATION
WAVELENGTH GIVES THE
SURGEON FULL FLEXIBILITY
WHEN CARRYING OUT INTRICATE
AIRWAY PROCEDURES.





The most feature-rich and highperformance airway scope in our collection.

Suitable for road or clinic-based vets, looking for superb image quality, durability and performance.



Unique infra-red laser filter - No white out



Titanium distal tip for durability



6 o'clock channel for midline surgical laser resection



German engineered

### DV3 LASER STAND

FOR STANDING
UPPER AIRWAY SURGERY

THE UNIQUE DESIGN OF THE SCOPE HOLDER GIVES GREAT STABILITY AND ALLOWS ME TO CONCENTRATE ON DEFLECTING THE SCOPE AND CUTTING TISSUE WITH THE LASER.



- ✓ Height adjustable for all breeds
- ✓ Optimal stability
- Lightweight 2KG







#### FOR CLINIC/HOSPITALS

The premium quality of the MEDView DV3 Airway combined with the MEDView PRO Endoscopy Trolley allows you to perform complex airway investigations, treatments and surgery in the clinic or hospital.

The MEDView PRO Endoscopy Trolley with 27-inch touchscreen enables you to work efficiently whether you have one scope or the full range.

Whether you're a university or a clinic-based hospital, the complete flexibility is there to perform airway, gastroscopy and oral dental examinations from one endoscopy station.

The **MEDView HDMI Pro Box** with MEDVos software acts as the link between your scope and screen. The large screen enables you to show split screen comparison to veterinarians and clients.

The trolley can store the full complement of MEDView endoscopes, instruments and accessories.







Real-time split screen





#### **MEDView DV3 Airway - Technical Specification**


Diameter	9mm	
Working length	1.6m	
Working channel	2.8mm	
Battery life	6 hours, mains power back up supply	
Portable battery accessory	Battery insufflator for uterine work	

Angulation	4 way   U: 140 D:110 L:110 R:110
Weight	9KG (including case)
Software	MEDVos
File format	JPEG/MPEG4
WIFI & Bluetooth	Built in

# Lemnardo DUAL 45



www.medequus.com

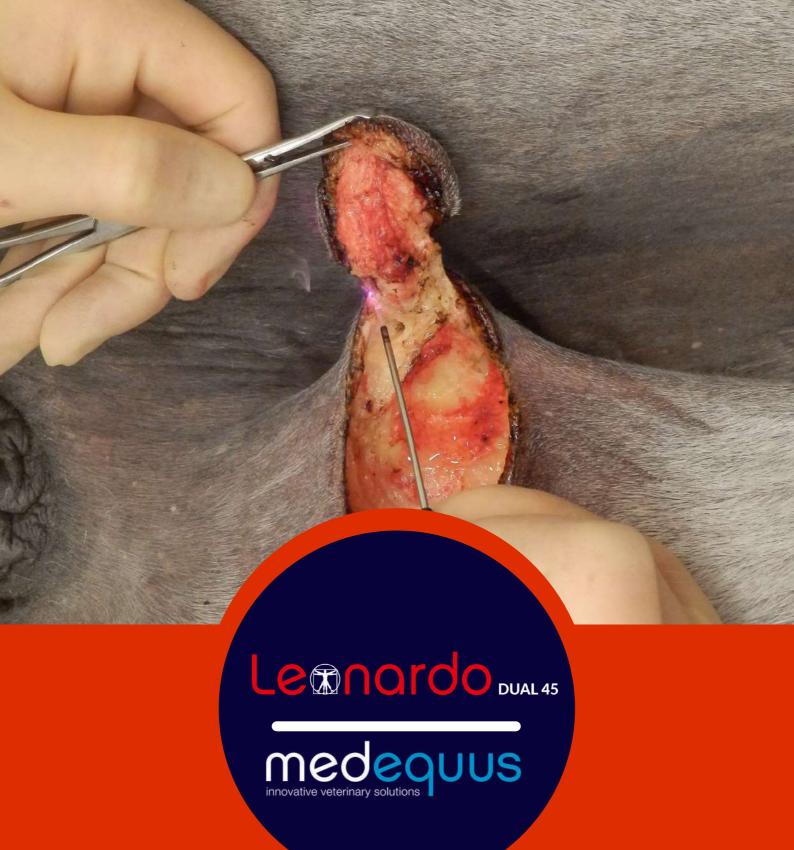

EU: +44 (0)1691 662 527 US: +1 (470) 5081969

**™** info@medequus.com

www.medequus.com

